

OPERATING
ROOM
ESSENTIALS

THE COMPLETE
SOLUTION
YOU CAN TRUST



VersaOne™
Bladed Trocar
VersaOne™ Access System

Medtronic
Further. Together

**EFFICACY
& SIMPLICITY.
EVERY TIME.**

The VersaOne™ Bladed Trocar is available in 5 mm, 11 mm, and 12 mm sizes

Smoother Insertion

Unique "dolphin nose" tip is designed to effectively dilate abdominal tissue layers with minimal insertion force

Port Fixation and Security

Advanced fixation ribs keep cannula from moving during surgical procedure


Protected Insertion

Parabolic shield is designed to retract over the blade once in the peritoneal cavity



Safety Lock

indicates when the blade is shielded and exposed

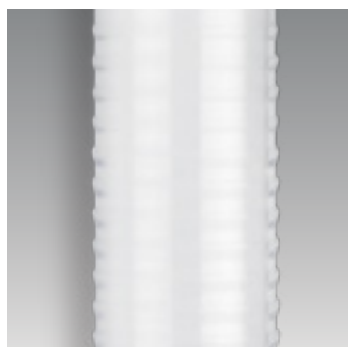

Proprietary ribbed cannula

for better port fixation and security compared to Applied Kii™ 12 mm trocar¹



Uniquely designed blade

divides tissue cleanly and precisely with control upon insertion through the abdominal wall



Enhanced seal system

allows for easy instrument exchange²

VersaOne™ Bladed Trocar

Sold six units per box

	5 mm Smooth	5 mm Fixation	11 mm Smooth	11 mm Fixation	12 mm Smooth	12 mm Fixation
Short Length Cannula (70 mm)	B5SHS	B5SHF	–	–	–	–
Standard Length Cannula (100 mm)	B5STS	B5STF	B11STS	B11STF	B12STS	B12STF
Long Length Cannula (150 mm)	–	–	–	–	B12LGS	B12LGF

VersaOne™ Universal Cannula

Sold six units per box

	5 mm Fixation	11 mm Fixation	12 mm Smooth	12 mm Fixation
Short Length Cannula (70 mm)	UNVCA5SHF	–	–	–
Standard Length Cannula (100 mm)	UNVCA5STF	UNVCA11STF	UNVCA12STS	UNVCA12STF
Long Length Cannula (150 mm)	UNVCA5LGF	–	–	UNVCA12LGF

Contact your Medtronic representative today for more information about the VersaOne™ Access System or visit medtronic.com/covidien.

1. Based on internal test report #2143-114. Applies to 12 mm trocars, when compared to Kii™* (12 mm, Z-thread cannula). March 2013.

2. Based on internal test report #2143-123. Applies to Medtronic 12 mm trocars, when compared to the Versaport™ Bladeless 12 mm trocar. March 2013.

© 2016 Medtronic. All rights reserved. Medtronic, Medtronic logo and Further, Together are trademarks of Medtronic.

™* Third party brands are trademarks of their respective owners. All other brands are trademarks of a Medtronic company.

16-eu-versaone-bladed-sell-sheet-906100

Medtronic



To contact us, please visit
medtronic.com/covidien/support/emea-customer-service

Use scan app to read