

CARDIOINSIGHT™ MAPPING VEST

Product specifications

Ordering Information

Order Number

CITVST0001S	Vest Size 1
CITVST0002M	Vest Size 2
CITVST0003L	Vest Size 3
CITVST0004XL	Vest Size 4

Vest specifications


Recommended Size

The vest comes in four sizes: 1, 2, 3, and 4.
Refer to Size Selection Reference Guide for information on vest sizing.

Material

Vest Panel	Cotton/spandex
252 Electrodes	Silver/Silver Chloride (AG/AGCl)
Gel	7% K-Cl Conductive Gel
Insulation Layer	Polyester
Connectors	Polycarbonate/ABS blend

Vest Cross-section

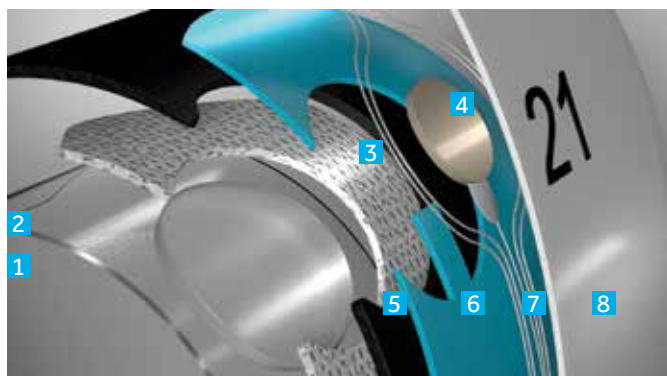

1. Liner 2. Conductive Gel 3. Adhesive 4. Sensor 5. Spandex 6. Dielectric
7. Circuit Trace 8. Polyester



The CardioInsight™ Noninvasive 3D Mapping System is intended for gathering, analysis, display, and storage of cardiac electrophysiological data and maps for analysis by a physician.

The CardioInsight™ Mapping vest is used with the CardioInsight™ Workstation. Together, they function as a noninvasive mapping system that collects chest ECG signals and combines these signals with CT scan data to produce and display simultaneous bi-atrial and biventricular, three dimensional cardiac maps.

The vest includes:

- 3 panels: front left, front right, and back
- 252 electrodes
- Connectors: connects the vest to the system in order to record cardiac activity
- The only materials in contact with the skin are the spandex, adhesive layer and electrode gel

The CardioInsight™ Mapping Vest is supplied non-sterile and packaged individually. It is a single use device.

Medtronic

CardioInsight™ Noninvasive 3D Mapping System



Essential Product Information

Indications and relevant warnings and instructions for use can be found in the product labeling supplied with each device.

Brief Statement

See the device manual for detailed information regarding the instructions for use, indications, contraindications, warnings, precautions, and potential adverse events. For further information, contact your local Medtronic representative or consult the Medtronic website at www.medtronic.com.

Medtronic

Europe

Medtronic International Trading Sàrl.
Route du Molliou 31
Case postale
CH-1131 Tolochenaz
www.medtronic.eu
Tel: +41 (0)21 802 70 00
Fax: +41 (0)21 802 79 00

United Kingdom/Ireland

Medtronic Limited
Building 9
Croxley Park
Hatters Lane
Watford
Herts WD18 8WW
www.medtronic.co.uk
Tel: +44 (0)1923 212213
Fax: +44 (0)1923 241004

UC201604043aEE © Medtronic 2017.
All rights reserved. Printed in Europe.