


Medtronic

FlexCath Cross™ Transseptal Solution

The FlexCath Cross transseptal solution combines a vessel dilator and transseptal needle into a single system. The FlexCath Cross system consists of an elongated shaft with a tapered tip and hollow stainless-steel needle to track over a guidewire in a similar fashion to conventional vessel dilators. The needle and shaft are connected to a handle incorporating a slide button to facilitate transseptal puncture. The needle lumen will accept guidewire sizes up to 0.032 inches in diameter.

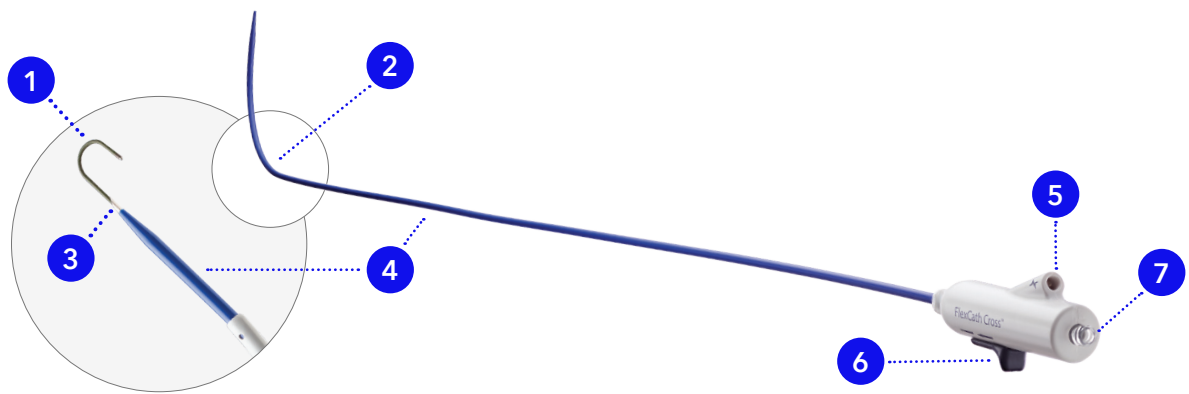

Contents

FlexCath Cross transseptal solution configurations include:

- FlexCath Cross Integrated Dilator/Needle
- Guidewire: 0.032 in, 180 cm, 3 mm J-tip, superstiff
- Electrosurgery and EGM Adapter Cables

Material

High-density polyethylene

Parts

1. Guidewire
2. Conforms to any sheath curve
3. Integrated needle (advances 4 mm)
4. Dilator shaft
5. EGM or RF energy connector
6. Needle slide button
7. Luer fitting guidewire lumen


Compatibility

- Guidewire compatibility: Up to 0.032 inch
- The distal tip of the handle is configured to mate with specific compatible transseptal vascular introducers, providing secure attachment during use. The FlexCath Cross system is intended to facilitate transseptal puncture when delivered through the compatible introducers in the table shown on the next page.




Ordering information and specifications

FlexCath Cross™ Transseptal Solution

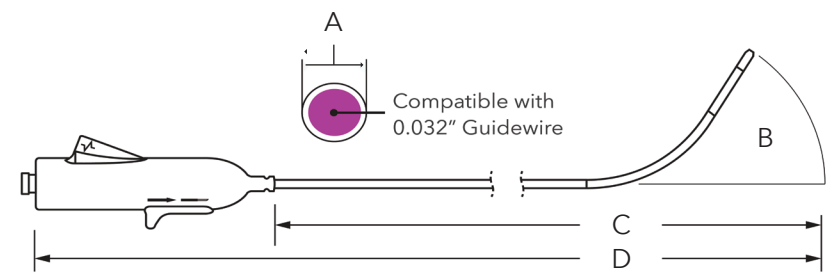

Product name	Compatible sheath	Order number	A - Outer diameter	B - Curve angle	C - Usable length	D - Total device length
FlexCath Cross – AG 61 cm	Agilis™ NxT – 61 cm	900300	8.5 Fr	70°	83.5 cm	92.4 cm
FlexCath Cross – AG 71 cm	Agilis NxT – 71 cm	900301	8.5 Fr	70°	93.4 cm	102.3 cm
FlexCath Cross – SL 63 cm	Swartz™ SL1 – 63 cm	900302	8.5 Fr	70°	67.7 cm	76.6 cm
FlexCath Cross – SL 81 cm	Swartz SL1 – 81 cm	900309	8.5 Fr	70°	85.6 cm	94.5 cm
FlexCath Cross – VZ 71 cm	CARTO VIZIGO® – 71 cm	900303	8.5 Fr	70°	93.4 cm	102.3 cm
FlexCath Cross – FC 65 cm	FlexCath Advance™ Sheath – 65 cm	900304	12 Fr	70°	83.5 cm	92.4 cm
FlexCath Cross – MH 63 cm	Merit Medical ML1 HeartSpan® Sheath – 63 cm	900306	8.5 Fr	70°	67.3 cm	76.2 cm
FlexCath Cross – MH 81 cm	Merit Medical ML1 HeartSpan Sheath – 81 cm	900307	8.5 Fr	70°	85.5 cm	94.4 cm
FlexCath Cross – WT 75 cm	WATCHMAN™ TruSeal™/ FXD Curve™ Access System – 75 cm	900308	12 Fr	70°	84.5 cm	93.4 cm

Mini Sheath, 65cm - L2 with FlexCath Cross™

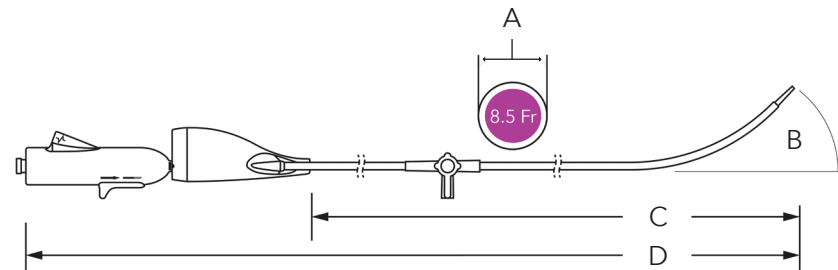

Product name	Compatible sheath	Order number	A - Outer diameter	B - Curve angle	C - Usable length	D - Total device length
Mini Sheath with FlexCath Cross Transseptal Solution	65 cm, 90 degrees, 8.5 Fr ID, 11 Fr OD	012852	11 Fr	90°	65 cm	84 cm

Medtronic

99 Hereford Street
Brampton, Ontario, L6Y 0R3
Toll-free: 800.268.5346
Tel: 905.460.3800

medtronic.ca

For a listing of indications, contraindications, precautions, warnings, and potential adverse events, please refer to the Instructions for Use.

©2024 Medtronic. All rights reserved. Medtronic, Medtronic logo and Engineering the extraordinary are trademarks of Medtronic. TM*Third party brands are trademarks of their respective owners. All other brands are trademarks of a Medtronic company. UC202311106bEC CA-CTL-0306-E Rev. 03/2024