


**Medtronic**

# Penditure™ LAA Exclusion System



# The Medtronic Penditure left atrial appendage (LAA) clip

- Fabric-free design for reduced risk of inflammation<sup>1</sup>
- True tip-first closure designed to prevent tissue from expressing beyond the distal tip of the clip<sup>2</sup>
- Recapturable, redeployable, and repositionable after deployment to ensure precise placement<sup>2</sup>
- Wide opening at the base offers good visibility, which may assist in proper clip placement

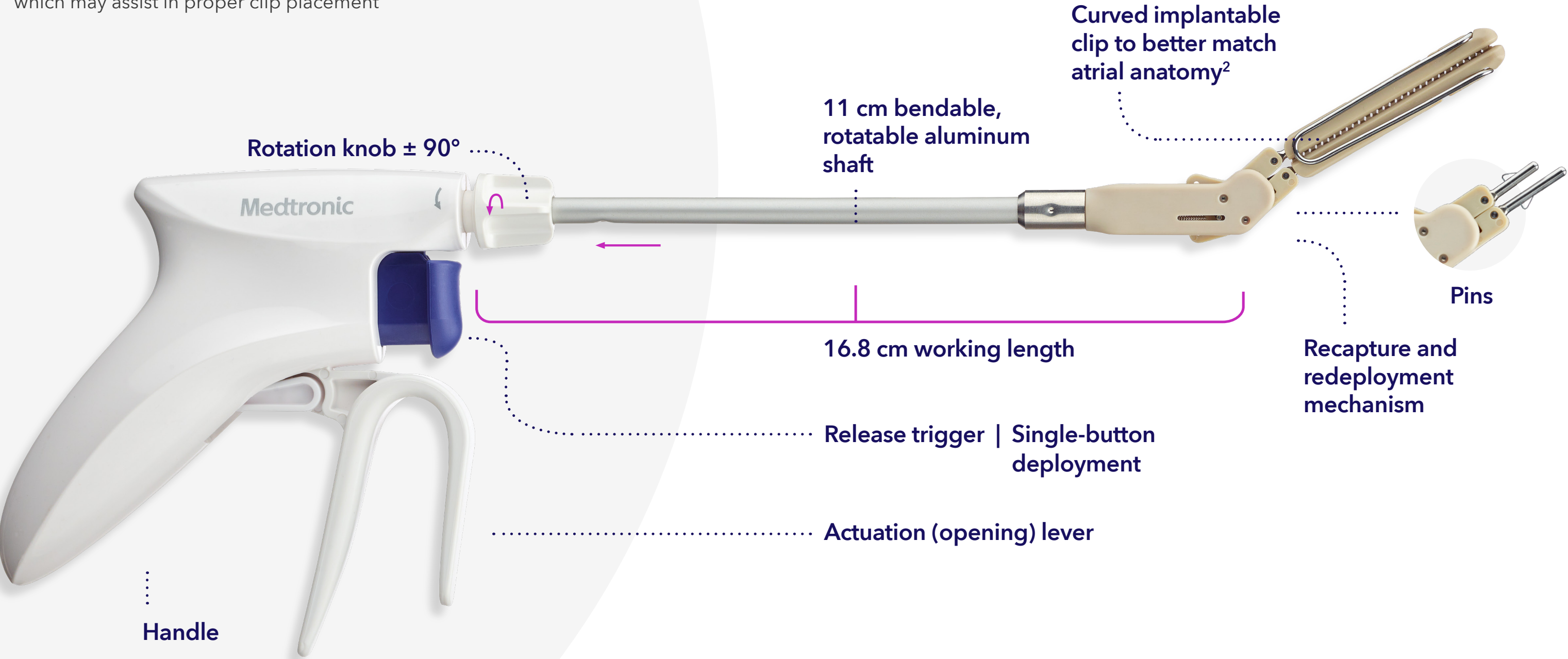


# Ordering information | Penditure LAA exclusion system

CFN	Length (A)	Height (B)	Width (C)	Opening (D) <sup>†</sup>
LAAC35 35 mm	39.6 mm	7.1 mm	11 mm	19.0 mm
LAAC40 40 mm	44.6 mm	7.3 mm	11 mm	19.9 mm
LAAC45 45 mm	50.0 mm	7.9 mm	11 mm	23.2 mm
LAAC50 50 mm	55.0 mm	8.6 mm	11 mm	20.0 mm
LAACSIZER	Packaged separately (3/pk)			

<sup>†</sup> Implant tip opening is a function of the moment arm and force of the spring. The location where spring force acts is further from the delivery system with longer clips. The result is that tip openings vary, and may not be consistent, with clip size.



Penditure clips are textured PEEK material with a subtle curve (radius). Nitinol springs provide the constant closing force.

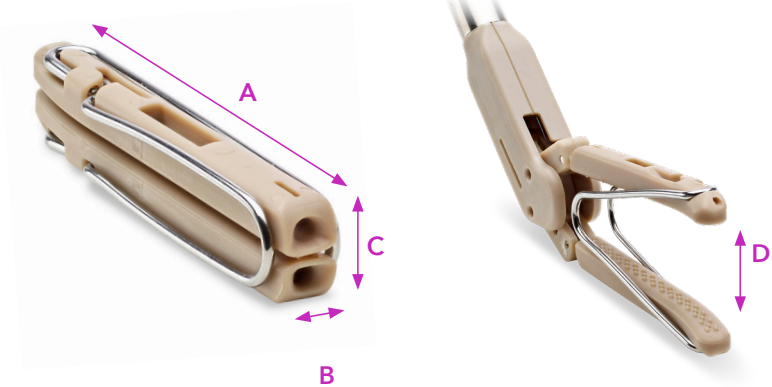

## Medtronic

99 Hereford Street  
Brampton, Ontario, L6Y 0R3  
Toll-free: 800.268.5346  
Tel: 905.460.3800

[medtronic.ca](https://www.medtronic.ca)

1. Medtronic animal testing on file, may not be indicative of clinical performance.
2. Medtronic data on file.

For a listing of indications, contraindications, precautions, warnings, and potential adverse events, please refer to the Instructions for Use.

©2025 Medtronic. All rights reserved. Medtronic, Medtronic logo and Engineering the extraordinary are trademarks of Medtronic. <sup>™</sup>\*Third party brands are trademarks of their respective owners. All other brands are trademarks of a Medtronic company. UC202306155bEC CA-CTL-0519-E Rev. 04/2025

