


Medtronic

Penditure™ LAA Exclusion System

Product matrix

Penditure LAA exclusion system consists of a single-use, sterile, self-closing, implantable clip preloaded on a single-use delivery system.

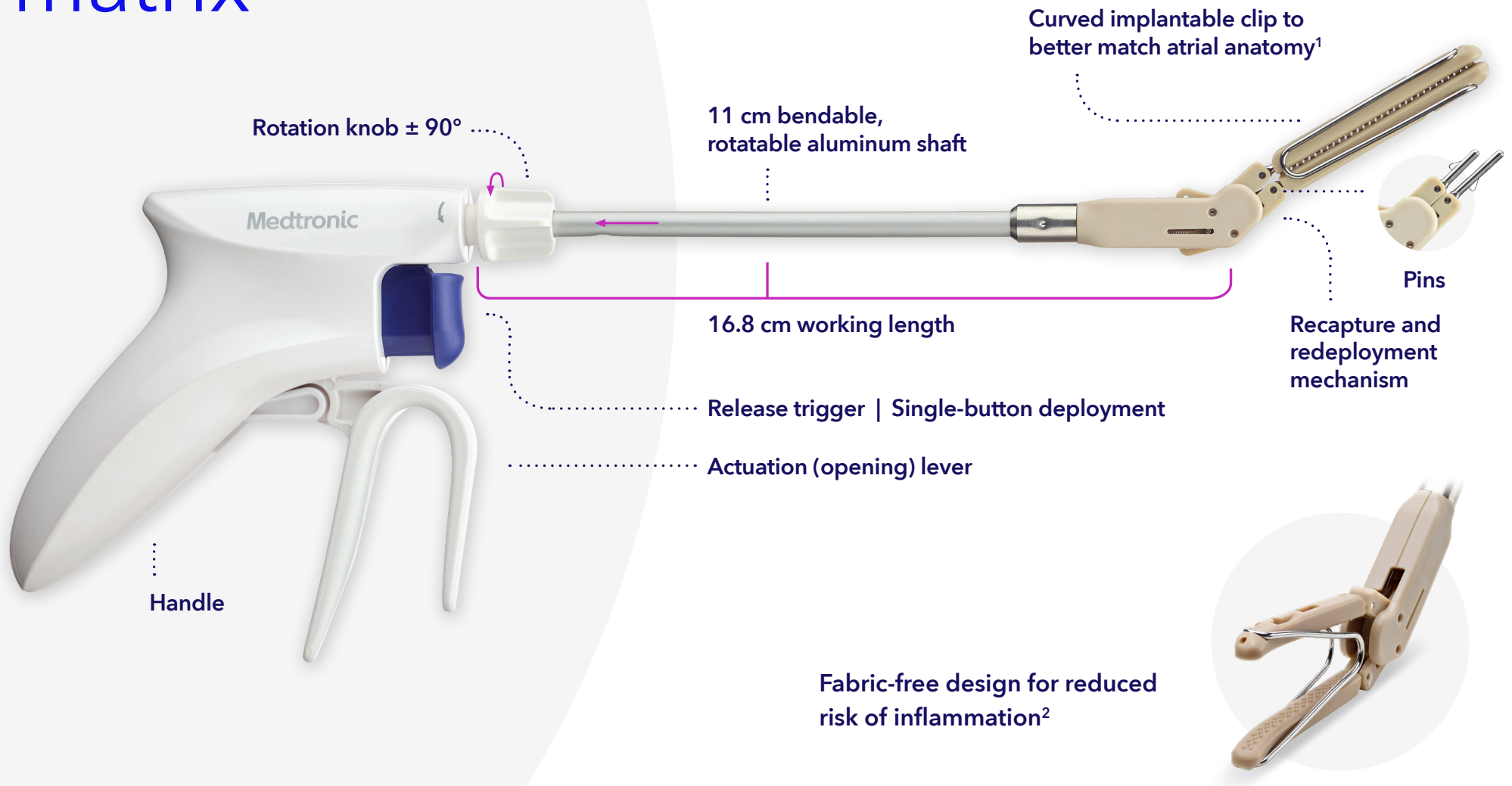

Penditure clip has a wide opening at the tip and at the base, designed to easily cross the appendage.¹

Clip size	Penditure clip tip opening when fully actuated†,‡	AtriClip™ Flex-V™ tip opening when fully actuated‡,§
35 mm	18.96 mm	Results pending
40 mm	19.89 mm	
45 mm	23.24 mm	
50 mm	19.95 mm	


Penditure clip fully actuated

A wide opening at the base offers good visibility, which may assist in proper clip placement.





† Implant tip opening is a function of the moment arm and force of the spring. The location where the spring force acts is further from the delivery system with longer clips. The result is that tip openings vary and may not be consistent with clip size.
‡ Design Characterization Report DO1035947.
§ Licensed as AtriClip FlexV.

Penditure clip fully actuated

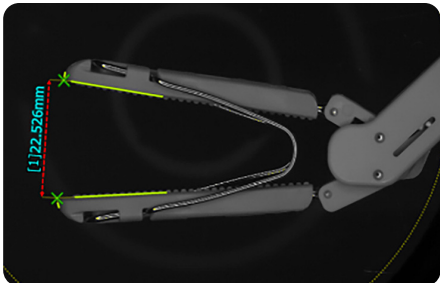


AtriClip Flex-V fully actuated§,1

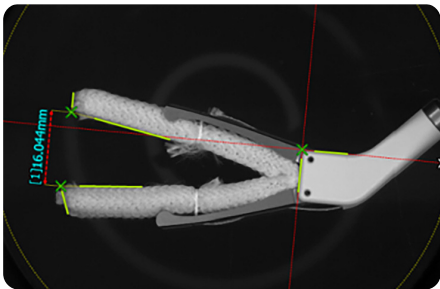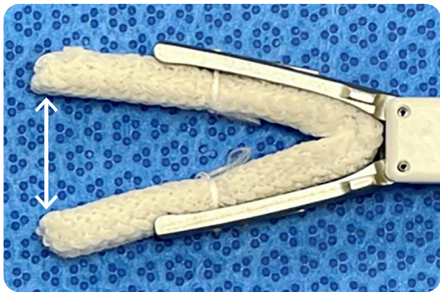


Wide u-shaped opening offers a large internal area.

Clip size	Penditure clip internal area when fully actuated [†]	AtriClip Flex-V internal area when fully actuated ^{†,‡}
35 mm	552.6 mm ²	Results pending
40 mm	523.6 mm ²	387.1 mm ²
45 mm	700.5 mm ²	406.1 mm ²
50 mm	729.6 mm ²	506.3 mm ²

Penditure clip fully actuated

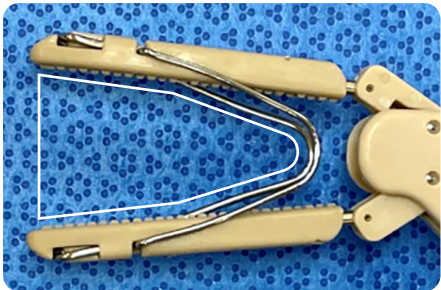


AtriClip Flex-V fully actuated

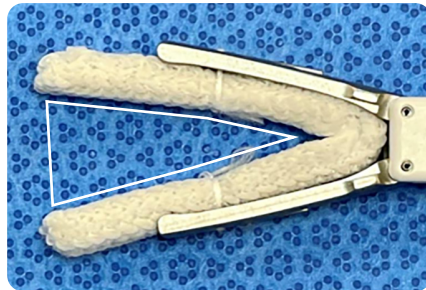

[†] Design Characterization Report DO1035947.

[‡] Licensed as AtriClip FlexV.


Penditure clip fully actuated

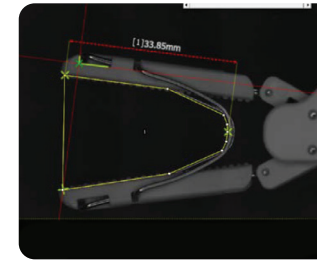



35 mm

AtriClip Flex-V fully actuated^{†,‡}



Results pending

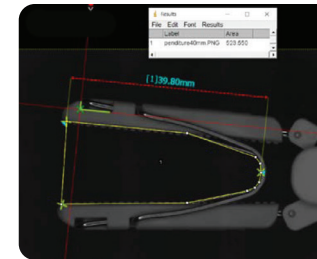


40 mm

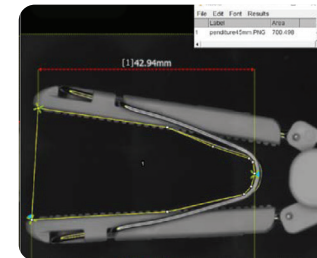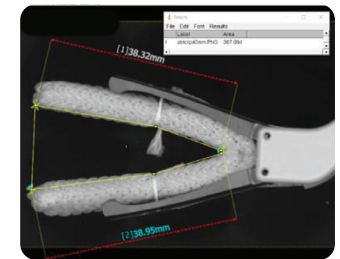

45 mm

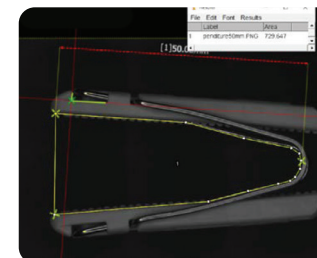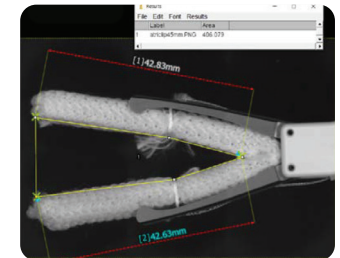

50 mm

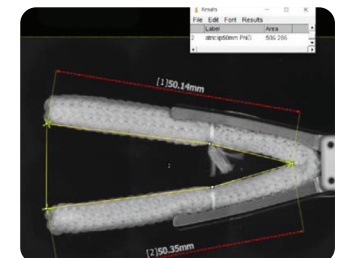

Tip-first closure helps provide more complete and secure exclusion of the LAA.³

Penditure LAA exclusion system was designed with true tip-first closure to prevent tissue from expressing beyond the distal tip of the clip.¹

Clip size	Penditure clip root opening [†]	AtriClip Flex-V root opening ^{†,‡}
35 mm	5.00 mm	Results pending
40 mm	5.27 mm	0 mm
45 mm	4.72 mm	0 mm
50 mm	5.22 mm	0 mm


Penditure clip
tip-first closure




AtriClip Flex-V
tip-first closure^{‡,1}



Penditure clip root opening

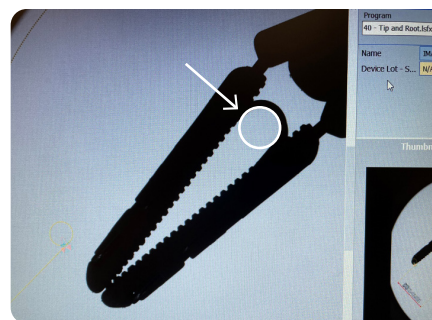

AtriClip Flex-V clip root opening^{‡,1}: N/A

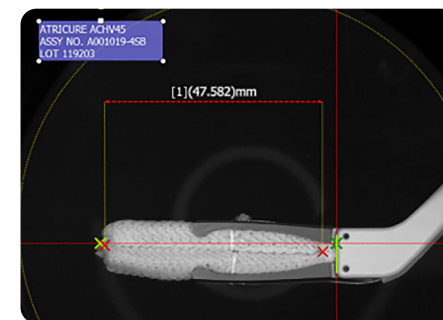



[†] Design Characterization Report DO1035947.

[‡] Licensed as AtriClip FlexV.

Product ordering information

CFN	Length (A)	Width (B)	Height (C)
LAAC35 35 mm	39.6 mm	11 mm	7.1 mm
LAAC40 40 mm	44.6 mm	11 mm	7.3 mm
LAAC45 45 mm	50.0 mm	11 mm	7.9 mm
LAAC50 50 mm	55.0 mm	11 mm	8.6 mm
LAACSIZER	Packaged separately (3/pk)		

Penditure clips are textured PEEK material with a subtle curve (radius). Nitinol springs provide the constant closing force.

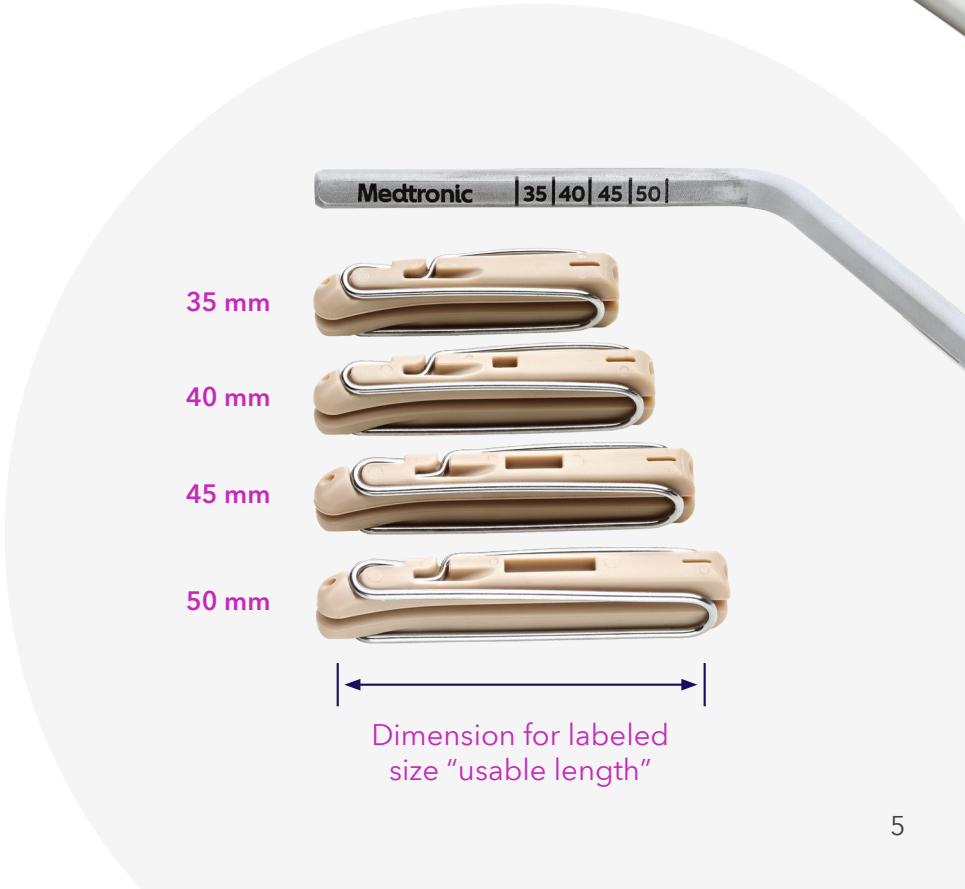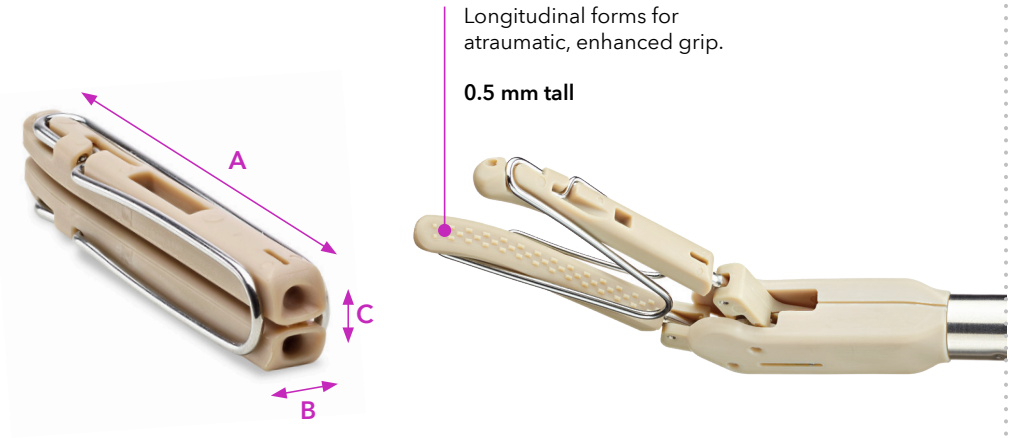

For a listing of indications, contraindications, precautions, warnings, and potential adverse events, please refer to the Instructions for Use.

References

1. Medtronic data on file.
2. Medtronic animal testing on file. May not be indicative of clinical performance.
3. Medtronic data on file. These tests may not be indicative of clinical performance.

Medtronic

99 Hereford Street
Brampton, Ontario, L6Y 0R3
Toll-free: 800.268.5346
Tel: 905.460.3800

[medtronic.ca](https://www.medtronic.ca)

©2025 Medtronic. All rights reserved. Medtronic, Medtronic logo and Engineering the extraordinary are trademarks of Medtronic. TM*Third party brands are trademarks of their respective owners. All other brands are trademarks of a Medtronic company.

UC202306157bEC CA-CTL-0520-E Rev. 05/2025