

Medtronic

FlexCath Cross™ Transseptal Solution

Zero exchange workflow





Efficiency

Zero exchange procedure

Zero wire and needle exchanges

Integrated needle/dilator design removes the need for needle and guidewire exchanges.


Zero sheath exchanges

Dilators compatible with leading industry sheaths eliminate the need for separate transseptal and therapeutic sheaths.




Efficiency Workflow

Traditional transeptal access workflow

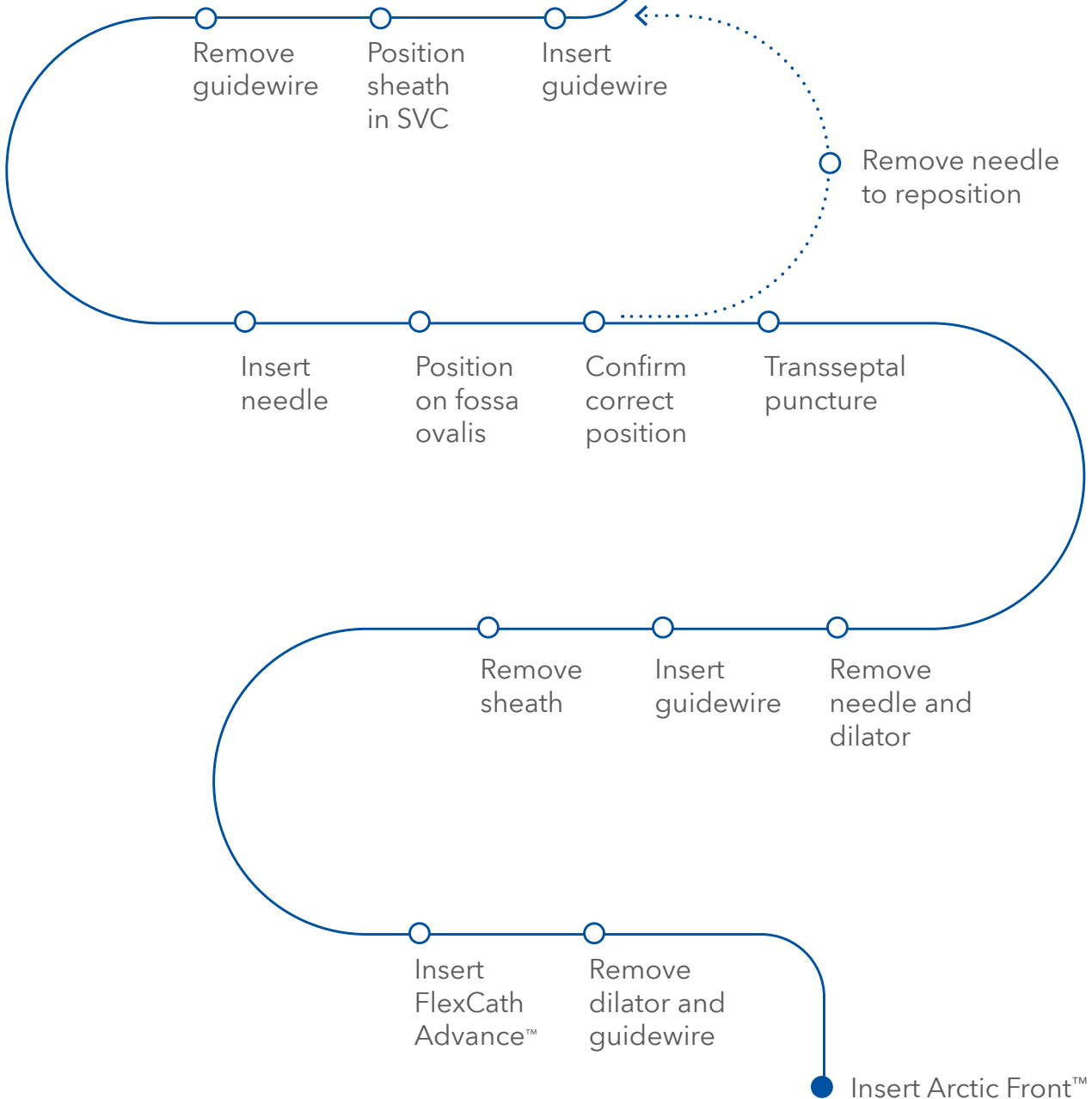

FlexCath Cross workflow with FlexCath Advance

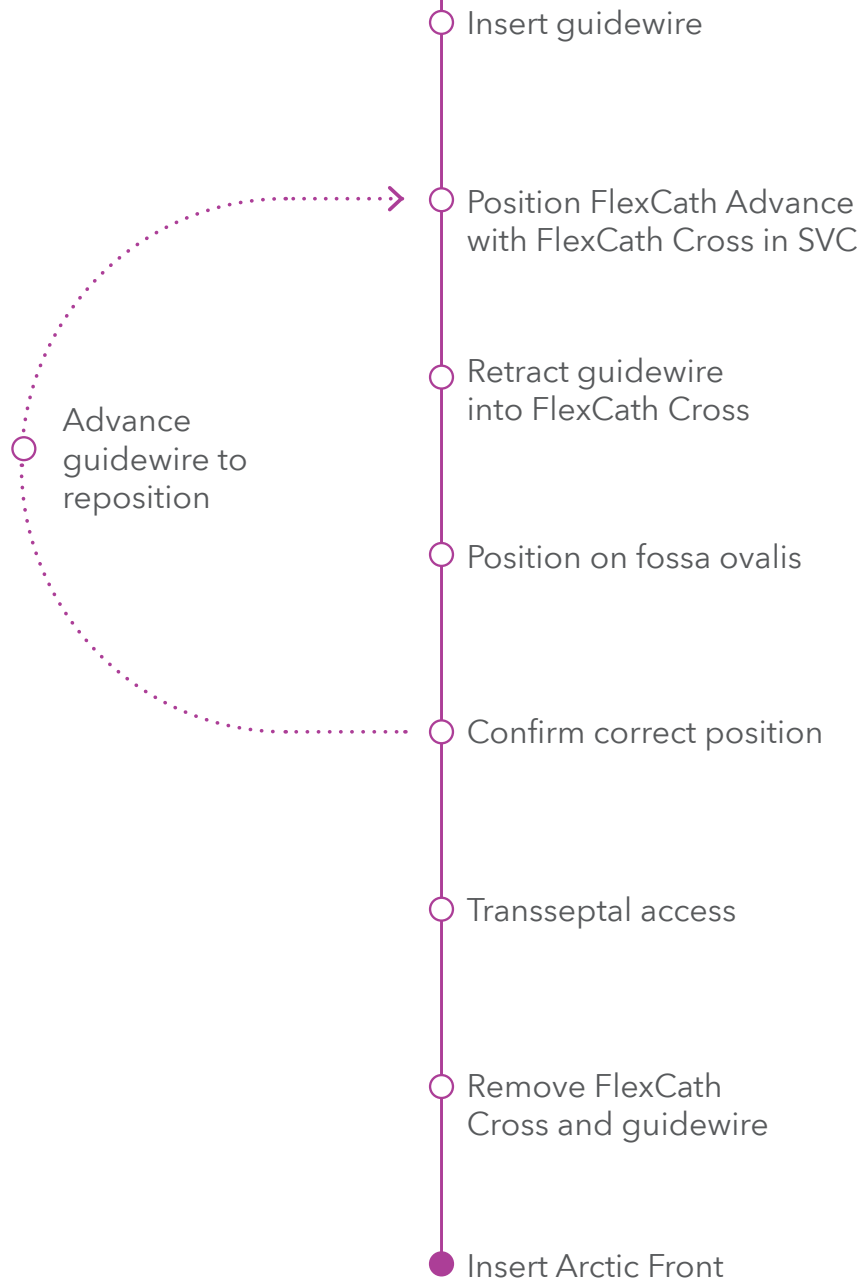

Versatile

Versatility for your procedures

Choose the FlexCath Cross model which fits with your existing sheath preferences for your left-heart procedures.

FlexCath Cross is the only transseptal access system on the market **cleared for both mechanical and RF crossing.**



Medtronic FlexCath Advance
Steerable Sheath

Compatible with industry-leading introducers



Abbott Agilis™ NxT
Steerable Introducer



Medtronic FlexCath Cross™
with Mini Sheath



Boston Scientific Watchman™
TruSeal™ and FXD™
Access Systems



Biosense Webster CARTO
VIZIGO®* Bi-directional
Guiding Sheath



Merit Medical HeartSpan®*
Sheath Introducer



Abbott Swartz™* Braided
Transseptal Guiding
Introducer, SL1

Controlled

Cross with control

The FlexCath Cross needle only deploys when activated, designed for ease of control during transseptal puncture.





Ordering info

Product name	Compatible sheath	Order number
FlexCath Cross – AG 61 cm	Agilis NxT – 61 cm	900300
FlexCath Cross – AG 71 cm	Agilis NxT – 71 cm	900301
FlexCath Cross – SL 63 cm	Swartz SL1 – 63 cm	900302
FlexCath Cross – SL 81 cm	Swartz SL1 – 81 cm	900309
FlexCath Cross – VZ 71 cm	CARTO VIZIGO – 71 cm	900303
FlexCath Cross – FC 65 cm	FlexCath Advance – 65 cm	900304
FlexCath Cross – MH 63 cm	Merit Medical HeartSpan Sheath – 63 cm	900306
FlexCath Cross – MH 81 cm	Merit Medical HeartSpan Sheath – 81 cm	900307
FlexCath Cross – WT 75 cm	Watchman TruSeal and FXD Access Systems – 75 cm	900308
FlexCath Cross with Mini Sheath – 65 cm	65 cm, 90 degrees, 8.5 Fr ID, 11 Fr OD	012852

Competitive matrix – transseptal needles

Feature	FlexCath Cross	BRK™*	SafeSept®*	NRG®*	VersaCross®*
Integrated needle/ dilator design reduces exchanges	✓	✗	✗	✗	✓
Spring-tensioned safety needle	✓	✗	✗	✗	✗
Indicated for use mechanically or with RF	✓	✗	✗	✗	✗
Reposition without requiring wire/needle exchange	✓	✗	✗	✗	✓
Engineered to pair with leading sheaths	✓	✗	✗	✗	✗

FlexCath Cross™ Transseptal Solution Indications (or Intended Use)

Prior to using these devices, please review the Instructions for Use for a complete listing of indications, contraindications, warnings, precautions, potential adverse events, and directions for use. FlexCath Cross devices are FDA cleared. Not all FlexCath Cross devices are CE marked.

Medtronic

710 Medtronic Parkway
Minneapolis, MN 55432-5604
USA

Toll-free in USA: 800.633.8766
Worldwide: +1.763.514.4000
medtronic.com

©2023 Medtronic. Medtronic, Medtronic logo, and Engineering the extraordinary are trademarks of Medtronic.™™Third-party brands are trademarks of their respective owners. All other brands are trademarks of a Medtronic company. UC202313327a EN 08/2023