

# RELIABILITY AND PERFORMANCE.

## MAHURKAR™\* Chronic Carbothane™\* Catheter




### PERFORMANCE

- 14.5Fr Double-D™ lumen design maximizes the cross-sectional area within the lumen
- Silicone extensions allow visibility of blood
- Ultem™ adapters withstand repeated use


### SAFETY

- Tal VenaTrac™ over-the-wire insertion stylets reduce risk of air embolism and blood loss<sup>1</sup>
- Safety tunneler allows for ease of insertion
- Large laser-printed volumes for easy reading

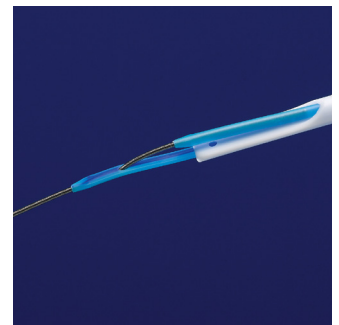

#### Tal VenaTrac™ stylets:

- Allow for sheathless catheter insertions<sup>1</sup>
- Promote over-the-wire catheter exchanges<sup>1</sup>

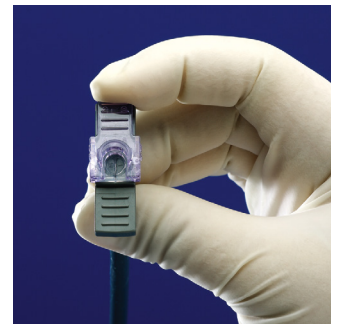

Valved pull-apart  
safety sheath

# ORDERING INFORMATION

## Priming Volumes

	Code	Description	Arterial (Red)	Venous (Blue)	Qty Case
<b>MAHURKAR™ Chronic Carbothane™ Catheter</b>					
SPORT PACK	8888119385	19/36 with sideholes	1.6 mL	1.7 mL	5
	8888123435	23/40 with sideholes	1.8 mL	1.9 mL	5
	8888128485	28/45 with sideholes	2.0 mL	2.1 mL	5
	8888133535	33/50 with sideholes	2.2 mL	2.3 mL	5
KITS	8888145251	19/36 with sideholes	1.6 mL	1.7 mL	5
	8888145252	23/40 with sideholes	1.8 mL	1.9 mL	5
	8888145253	28/45 with sideholes	2.0 mL	2.1 mL	5
	8888145254	33/50 with sideholes	2.2 mL	2.3 mL	5
<b>MAHURKAR™ Chronic Carbothane™ Catheter with Tal VenaTrac™ Over-the-Wire Insertion Stylets (One Pair)</b>					
SPORT PACK	8888119362	19/36 with sideholes	1.6 mL	1.7 mL	5
	8888123402	23/40 with sideholes	1.8 mL	1.9 mL	5
	8888128452	28/45 with sideholes	2.0 mL	2.1 mL	5
	8888133502	33/50 with sideholes	2.2 mL	2.3 mL	5
KITS	8888145027	19/36 with sideholes	1.6 mL	1.7 mL	5
	8888145028	23/40 with sideholes	1.8 mL	1.9 mL	5
	8888145029	28/45 with sideholes	2.0 mL	2.1 mL	5
	8888145030	33/50 with sideholes	2.2 mL	2.3 mL	5

Components	Qty	Sport pack	Kit
MAHURKAR™ Chronic Carbothane™ Catheter	1	■	■
16Fr valved pull-apart safety sheath	1	■	■
Tunneling stylet, 21 cm	1	■	■
12Fr tissue dilator	1	■	■
14Fr tissue dilator	1	■	■
Sealing caps	2	■	■
18 G x 7 cm introducer needle	1		■
J/Straight .038" x 120 cm guidewire	1		■
No. 11 scalpel	1		■
Telfa™ Island Dressings	2		■
4" x 4" cotton gauze sponges	4		■
12 mL syringe	1		■

For more information, visit us at [medtronic.com/covidien/products/dialysis-catheters](https://www.medtronic.com/covidien/products/dialysis-catheters)

1. Mojibian HR, Oneal P, Eliseo D, Tal MG. VenaTrac™ device for over-the-wire placement of chronic dialysis catheters. *J Vasc Access*. 2010;11(4):342-345.

#### MAHURKAR™ Chronic Carbothane Catheter

The Mahurkar 14.5 Fr Chronic Carbothane catheter is intended for acute and chronic hemodialysis, apheresis, and infusion. It may be inserted either percutaneously or by cutdown.

Do not use the catheter in thrombosed vessels or for subclavian puncture when ventilator is in use.

See the device manual for detailed information regarding the implant procedure, indications, contraindications, warnings, precautions, and potential complications/adverse events. For further information, please call Medtronic, in the United States call 1-800-962-9888 and for international inquiries, call 508-261-8000, and/or consult [www.medtronic.com/covidien](https://www.medtronic.com/covidien)

© 2020 Medtronic. All rights reserved. Medtronic, Medtronic logo and Further, Together are trademarks of Medtronic.™ Third party brands are trademarks of their respective owners. MAHURKAR™ is a U.S. registered trademark of Sakharam D. Mahurkar, used under license. All other brands are trademarks of a Medtronic company. 09/2020 – US-RC-2000026 – [WF#1906407]

15 Hampshire Street  
Mansfield, MA 02048

800.962.9888  
508.261.8000

[medtronic.com/covidien](https://www.medtronic.com/covidien)

# Medtronic