




The DEHP-free Shiley™ oral/nasal endotracheal tube  
with intermediate cuff

# For pediatric through adult patients.

The DEHP-free Shiley™ oral/nasal endotracheal tube with intermediate cuff can be used for pediatric through adult patients.

- 5.0-10.0 mm size options feature an intermediate high-volume, low-pressure cuff
- 3.0-4.5 mm size options feature an intermediate low-profile cuff to meet the unique needs of pediatric patients

Compared to the legacy version, this product is made with non-DEHP PVC material, which softens at body temperature and molds to the airway.

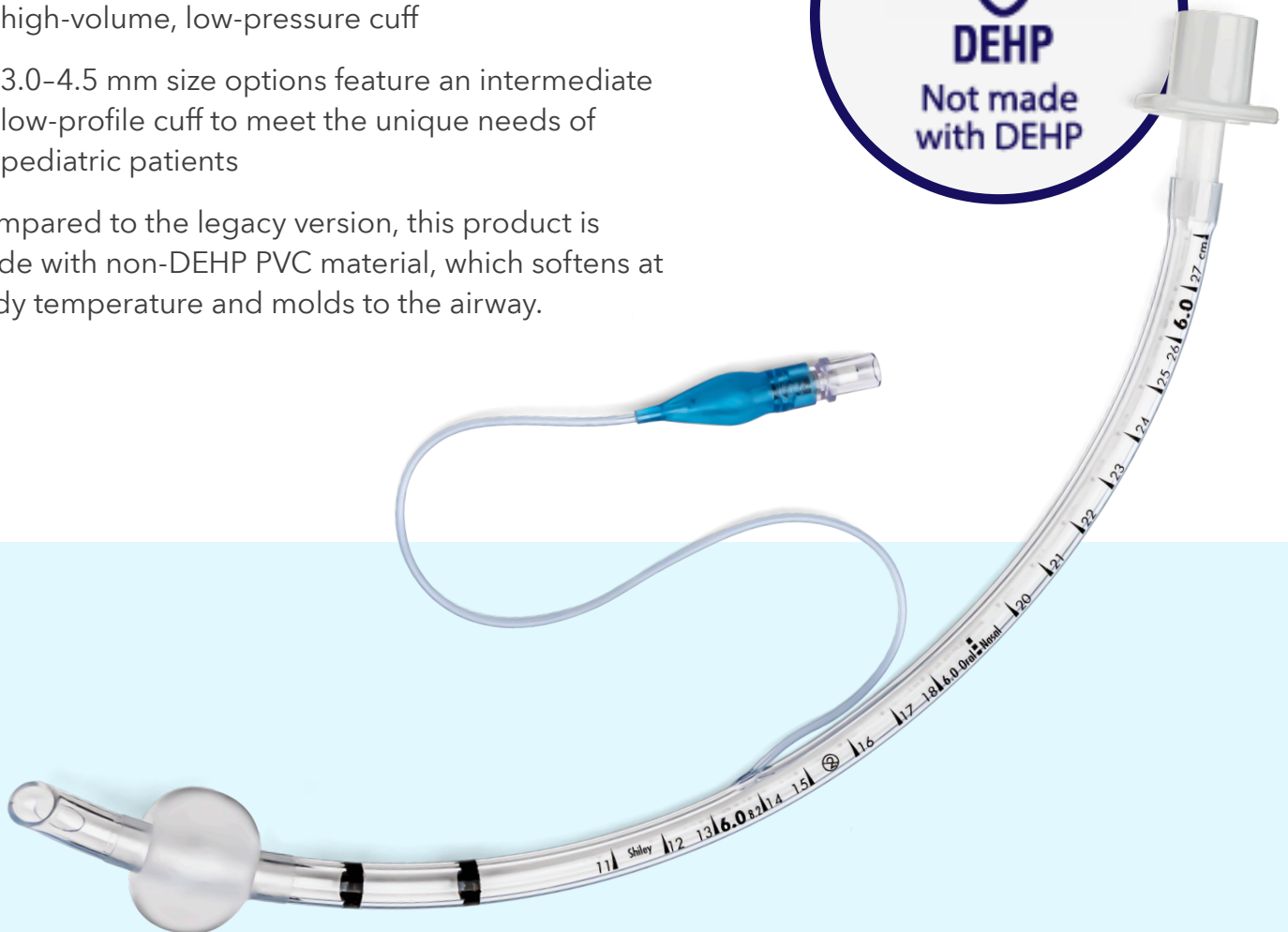

# Features and benefits

## Sizes 5.0-10.0

### A standard 15 mm connector

For connection to respiratory and anesthesia equipment

### Depth marks in centimeter increments

Printed along a portion of the length of the device

### Magill curve

For easier tube insertion

### Embedded radiopaque filament within its wall

Present down the length of the tube and visible on x-ray

### Glottic print marks – full circumference

For tube visualization during intubation and tube placement when the tip is out of view and beyond the vocal cords†

### Intermediate high-volume, low pressure cuff

Less material compared to the original Hi-Lo and Curity™ product

### Hooded tip

To reduce tip sharpness

### Beveled end with Murphy eye

To aid in visibility and insertion through the vocal cords and into the trachea

## Sizes 3.0-4.5

### Glottic print marks – 180°

### Intermediate low profile cuff

Elongated, low-profile cuff used in pediatrics

### Straight tip with a Murphy eye

## Ordering Information

| CFN   | I.D. (mm) | O.D. (mm) | Length (mm) | Cuff Ø (mm) |
|-------|-----------|-----------|-------------|-------------|
| 87430 | 3.0       | 4.4       | 160         | 7.8         |
| 87435 | 3.5       | 5.0       | 180         | 9.6         |
| 87440 | 4.0       | 5.7       | 200         | 9.6         |
| 87445 | 4.5       | 6.3       | 220         | 13.0        |
| 87450 | 5.0       | 6.9       | 240         | 20.9        |
| 87455 | 5.5       | 7.5       | 270         | 21.8        |
| 87460 | 6.0       | 8.2       | 280         | 23.9        |
| 87465 | 6.5       | 8.9       | 290         | 24.7        |
| 87470 | 7.0       | 9.5       | 300         | 28.2        |
| 87475 | 7.5       | 10.2      | 310         | 28.7        |
| 87480 | 8.0       | 10.8      | 320         | 32.3        |
| 87485 | 8.5       | 11.4      | 330         | 32.9        |
| 87490 | 9.0       | 12.1      | 335         | 34.0        |
| 87495 | 9.5       | 12.8      | 335         | 35.1        |
| 87410 | 10.0      | 13.5      | 335         | 35.2        |

Carton quantity: 10. All tubes are delivered sterile packed, for single use only.

† As indicated by ISO 5361:2016.



Scan the QR code to learn more or contact your local Medtronic representative.

©2025 Medtronic. Medtronic, Medtronic logo, and Engineering the extraordinary are trademarks of Medtronic. All other brands are trademarks of a Medtronic company. 06/2025 – US-RE-2200448 – [WF# 7477881]

**Medtronic**