


Epsila EV™ sternal tunneling tool

Model EAZ101

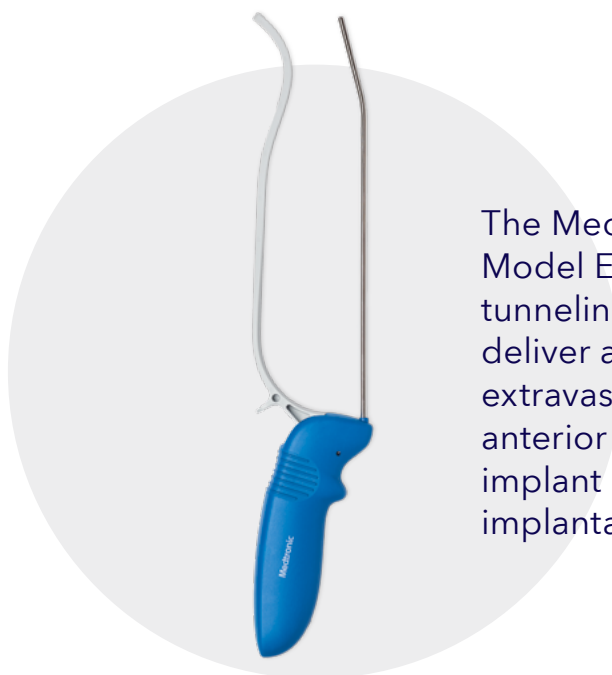


The Medtronic Epsila EV Model EAZ101 sternal tunneling tool is designed to deliver an introducer and an extravascular lead into the anterior mediastinum during implant of an extravascular implantable device system.

Product specifications

The sternal tunneling tool consists of the following components:

- 1 A **handle**.
- 2 A **thumb tab** that can be used to raise and lower the external guide.
- 3 A **stainless-steel tunneling rod** that delivers an introducer to the anterior mediastinum. The tunneling rod has a preformed bend and is malleable to accommodate patient anatomy.
- 4 An **external guide** that remains above the skin and indicates the distance and direction of the tunneling rod. The external guide is hinged and removable to accommodate physician preference and patient anatomy.

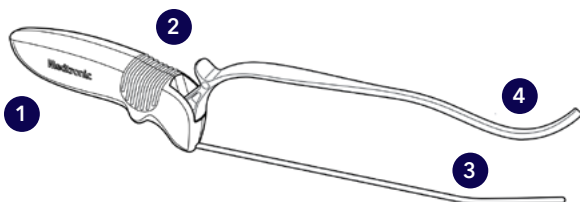

This tool is for single use only.

Physical characteristics

| | |
|------------------------|---------------------------------------|
| Length | |
| Overall length | 35.5 cm |
| Tunneling length | 22.3 cm |
| Tunneling rod diameter | 3.0 mm (9 Fr) |
| Material | |
| Handle | Polycarbonate |
| Tunneling rod | 304 stainless steel |
| External guide | Acrylonitrile butadiene styrene (ABS) |
| Recommended introducer | |
| Diameter | 3.0 mm (9 Fr) |
| Length | 19.13 cm |

Brief Statement

This material does not replace or supersede the instructions for use. It should not be considered the exclusive source of information, and should be used in conjunction with the device manual.

See the device manual for detailed information regarding the instructions for safe use, the implant procedure, indications, contraindications, warnings, precautions, and potential adverse events. If using an MRI SureScan™ device, see the MRI SureScan™ technical manual before performing an MRI. For further information, contact your local Medtronic representative and/or consult the Medtronic website at medtronic.eu.

For applicable products, consult instructions for use on medtronic.com/manuals. Manuals can be viewed using a current version of any major internet browser. For best results, use Adobe Acrobat® Reader with the browser.

Europe

Medtronic International Trading Sàrl.
Route du Molliau 31
Case postale
CH-1131 Tolochenaz
www.medtronic.eu
Tel: +41 0 21 802 70 00
Fax: +41 0 21 802 79 00

United Kingdom/Ireland

Medtronic Limited
Building 9
Croxley Park
Hatters Lane
Watford
Herts WD18 8WW
www.medtronic.co.uk
Tel: +44 0 1923 212213
Fax: +44 0 1923 241004

UC202103595 EE ©2023 Medtronic. All rights reserved. Medtronic, Medtronic logo, and Engineering the extraordinary are trademarks of Medtronic. Adobe and Acrobat Reader are registered trademarks of Adobe Systems incorporated in the United States and/or other countries. Printed in Europe.